### EASYGEL INX X100

Easier printing with gelatin

#### "Printing was easy and the user experience was great!" - Anni Mörö, COO @ STEMSIGHT Oy.

Faculty of Medicine and Health Technology, BIOMEDITECH, Tampere University, Finland




EASYGEL INX<sup>®</sup> X100 is a gelatin-based, shear thinning, cell-interactive, extracellular matrix (ECM)-mimicking ink. It provides all the benefits of gelatin with the addition of shear thinning behaviour, thereby allowing easy printing at 37 °C.

Based on gelatin derived from natural collagen, EASYGEL INX<sup>®</sup> X100 has been modified with photocrosslinkable functional groups meaning that it can be printed with unprecedented efficiency. It resembles the natural ECM and has exceptionally high cell viabilities.

#### **BIOLOGICAL APPLICATIONS**

EASYGEL INX<sup>®</sup> X100 has been used to generate and sustain 3D cellular structures of a variety of human tissues, including adipose derived stem cells (ASC), and human foreskin fibroblasts (HFF). For more information on the biological applications of EASYGEL INX<sup>®</sup> X100 and the parameters used to generate these 3D cellular structures, contact us on <u>info@bioinx.com</u>

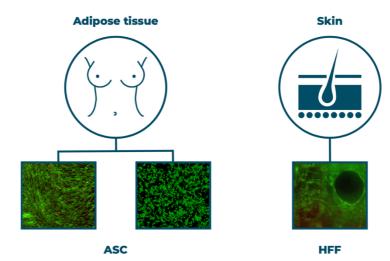

BIOINX

www.bioinx.com • info@bioinx.com

## EASYGELINX X100

#### Easier printing with gelatin

EASYGEL INX<sup>®</sup> X100 combines all the benefits of conventional gelatin / Gel-MA based inks with a highly improved printing process thanks to its shear thinning behavior. As a result, it can be printed in a straightforward and reproducible way at 37°C both in the presence or absence of cells (Figure 1). In this respect, it overcomes the narrow processing range limitations associated with printing a conventional Gel-MA based ink.



Conventional Gel-MA at 37°C vs. EASYGEL INX<sup>®</sup> at 37°C

Figure 1: Comparison of flow behavior of conventional Gel-MA based inks vs EASYGEL INX® X100 at 37°C

Since EASYGEL INX<sup>®</sup> X100 is gelatin-based, it is characterized by all the favorable gelatin properties including cell interactivity, biodegradability and the potential for cell encapsulation with high cell viability. It acts as an ideal mimic of the natural extracellular environment, making it suitable for a range of tissues in combination with unprecedented processing ease. After printing, the material can be photo-crosslinked, resulting in a physiologically stable hydrogel.

#### BENEFITS OF EASYGEL INX<sup>®</sup> X100

- ✓ Easy Printing Shear thinning behavior enables easy deposition at 37°C
- ✓ Biodegradability Enables cellular remodeling of the printed matrix
- ✓ Easy Handling Delivered in a ready-to-use cartridge
- ✓ UV-curable Efficient UV-based crosslinking
- ✓ Reproducibility Production under strict quality control

### EASYGELINX X100

Easier printing with gelatin

Table 1: Typical benefits of EASYGEL INX® over conventional bioinks

	Conventional Alginate-based bioink	Conventional gel-MA based bioink	EASYGEL
Printing at 37°C	(3)	(2)	$\bigotimes$
Cell-interactivity		$\bigotimes$	$\bigotimes$
Biodegradability		$\bigcirc$	$\oslash$
Biocompatibility	$\bigcirc$	$\bigcirc$	$\oslash$
Shear-thinning	$\bigcirc$	(2)	$\oslash$
Hydrogel		$\bigotimes$	$\oslash$

#### **PROPERTIES & PROCESSING**

EASYGEL INX<sup>®</sup> X100 is a transparent gel at room temperature. The physical characteristics of the product are listed in Table 2.

EASYGEL INX<sup>®</sup> X100 reveals a shear-thinning behavior as shown in Figure 2 (a). At high shear rates, it exhibits a low viscosity which allows easy injection from the printing nozzle. However, at low shear rates it has a high viscosity, which is required for shape retention after deposition. This is favorable for extrusion-based 3D printers as the ink should easily be injected through the printing nozzle and the post-injection flow should be minimized to prevent structural deformation. These properties are very rarely seen in conventional gelatin-based inks, making their processing rather complicated.

Physical Properties	EASYGEL INX <sup>©</sup> X100	
рН	6.5 – 8.5	
Viscosity (low shear) (Pa.s)	200 – 2500	
Viscosity (high shear) (Pa.s)	0.1 – 1	
Storage modulus after UV crosslinking (kPa)	2-6	

Table 2: Physical properties of EASYGEL INX<sup>®</sup> X100



### EASYGELINX X100

Easier printing with gelatin

To enable an optimal printing process ensuring shape fidelity, an ink should not only exhibit shear-dependent viscosity, but also it should regain its structure quickly upon deposition on the printing surface after being injected from the needle. The rate of structure recovery was studied via oscillatory step strain tests at strain values of 10%, 3000% and 10% in sequential order. As seen in Figure 2 (b), EASYGEL INX<sup>®</sup> X100 exhibits a rapid structure recovery when the external forces are reduced.

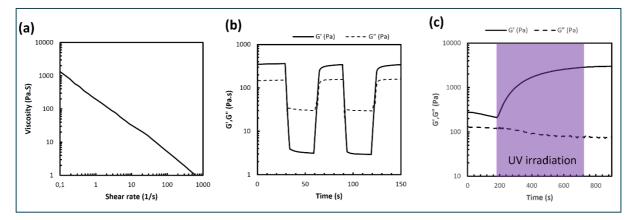


Figure 2: (a) Flow curve of EASYGEL INX® X100 as a function of shear rate, (b) Storage and loss modulus of EASYGEL INX® X100 tested via step strain tests and (c) storage and loss moduli of EASYGEL INX® X100 recorded via a rheometer during UV irradiation.

EASYGEL INX<sup>®</sup> X100 is photo-crosslinkable, and therefore the structures can be illuminated with UV light during or after the printing process. After switching on the UV light source, the ink exhibits a fast crosslinking reaction as indicated by the steep increase of storage modulus (Figure 2 (c)). At the end of the irradiation process, the ink reached a sufficient storage modulus for mechanical integrity while being ideal for soft tissue engineering.

#### **3D PRINTER COMPATIBILITY**

Our resins have been used repeatedly and successfully with the following printers:

- ✓ Regemat3D
- ✓ FelixBio
- ✓ Cellink BIOX

If you would like to discuss your printer's compatibility with our bioinks, please contact us at info@bioinx.com

